


Designed expressly for Viral Flu checkup

HiViral™ Transport Kit

Safe, Easy and Efficient

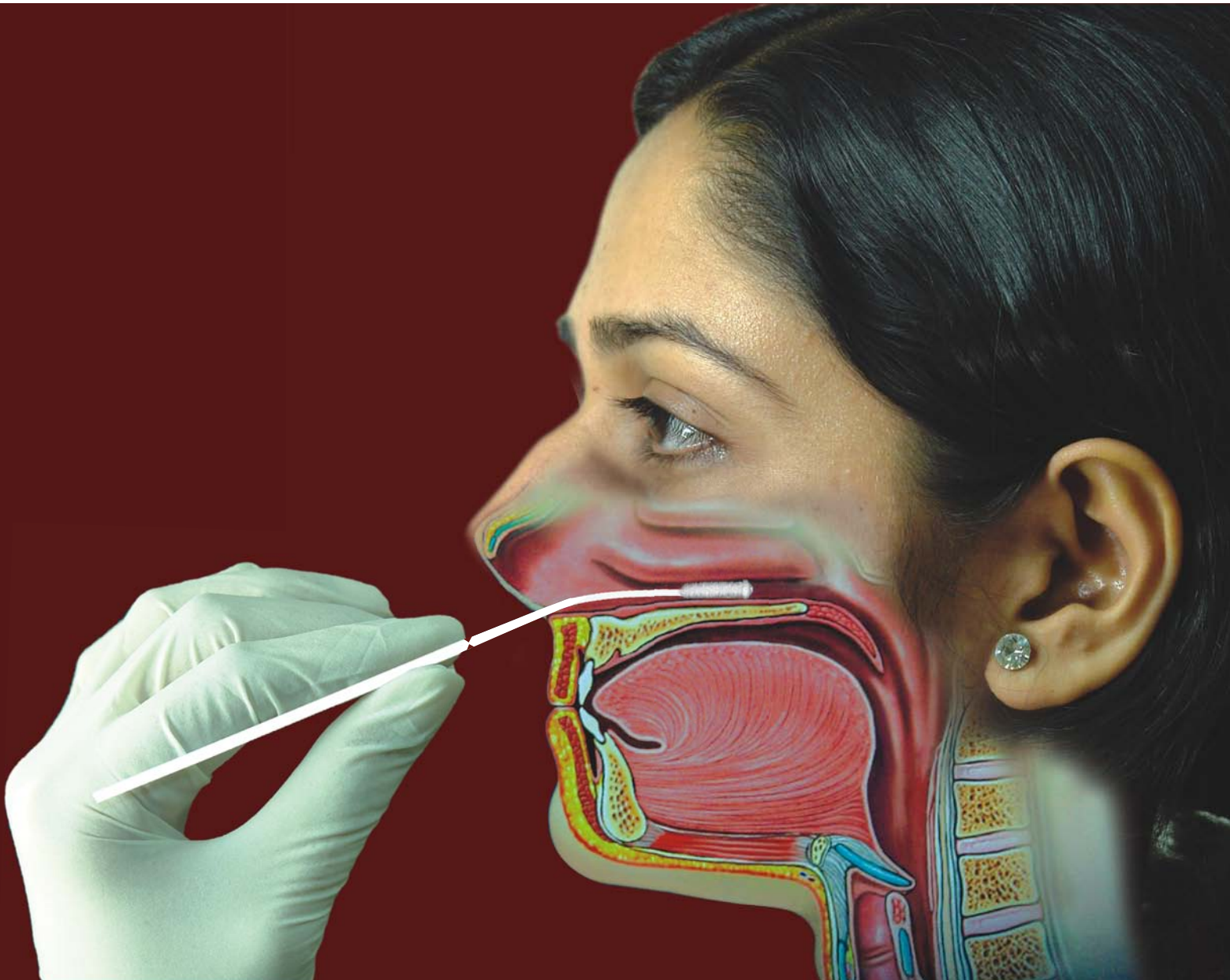

Designed expressly for Viral Flu checkup



HIMEDIA®


For life is precious

HiViral™ Transport Kit

MS1760C-50NO

Safe, Easy and Efficient

Combined system for Collection
and Transportation of Viruses

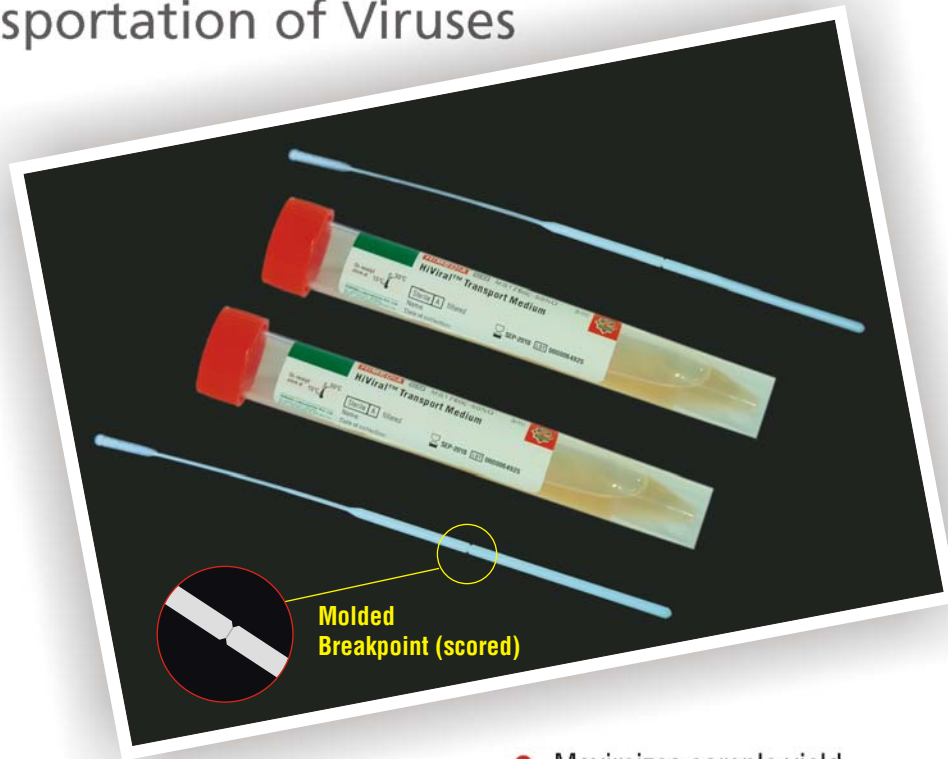

Self-contained Viral transport system,
designed to meet your viral testing needs

- Maximizes sample yield
- Maximizes recovery
- Compatible with cell culture and molecular based assays
- PW1172 has a molded breakpoint which allows swab to be broken into the tube

FEATURES in the Kit


HiViral™ Transport Medium

Designed to maintain viral viability and transport viruses in active form to the laboratory for isolation.

- ▶ Protein stabilizes the virus
- ▶ Buffers control the pH
- ▶ Antibacterial and Antifungal agents prevent bacterial and fungal contamination
- ▶ Phenol red acts as a pH indicator and ensures that medium is suitable for use at the time of specimen collection
- ▶ Cryoprotectant helps in preserving the viruses if samples are frozen for prolonged storage
- ▶ Dual temperature compatible, can be stored and transported at room as well as refrigerated temperature


Flocked Nylon Swab

- ▶ 6" long shaft enables easier sampling at difficult sites without direct exposure to the user
- ▶ Short perpendicular nylon fiber strands attached to flexible molded plastic results in efficient collection and release of particulate matter
- ▶ Recent references have shown better sample yield with flocked swabs than conventional swabs
- ▶ Higher sample yield maximizes sensitivity of cell culture and molecular techniques
- ▶ Flexible plastic shaft delivers better patient comfort
- ▶ Packed individually in a peel/tear to open pouch format

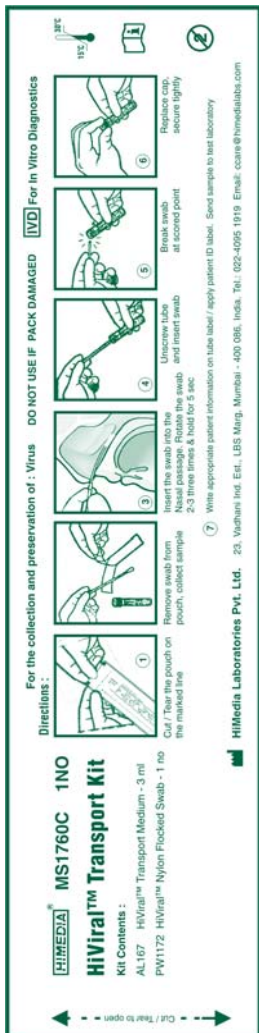


Maximum patient comfort with:-


I) HiMedia Flocked Nylon Swab for Nasopharyngeal Sample Collection

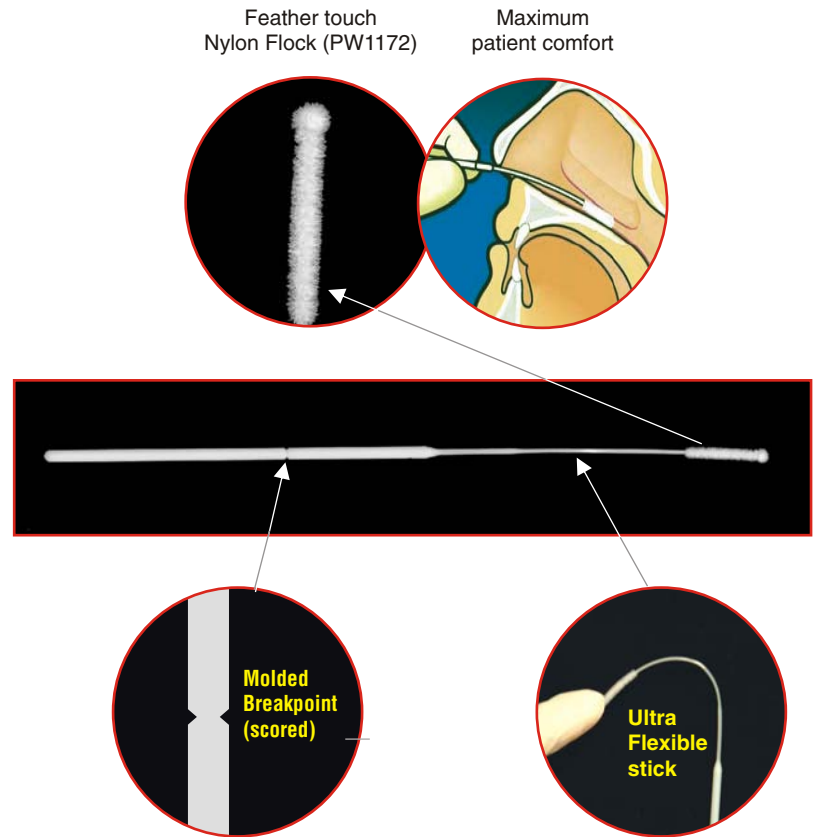

II) Throat Sample Collection

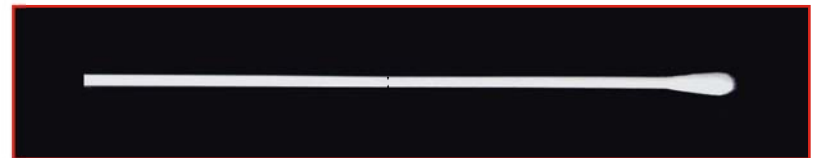


Throat (Polyester) scored Swab

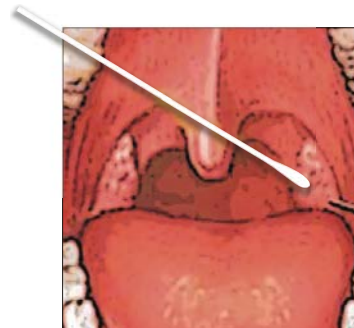


Suggested Throat sample collection area
(Sample may be collected from other areas of the throat as per user established protocols)

HIMEDIA[®]

For life is precious

Product Name	Code	
1. HiViral Transport Kit (C) [Nasopharyngeal swab + Medium (in self-standing tube)]	MS1760C	For MS1760C, MS1760F & MS1760G 
2. HiViral Transport Kit (F) [Throat (Polyester) scored swab + Medium (in self-standing tube)]	MS1760F	
3. HiViral Transport Kit (G) [Foam tipped swab + Medium (in self-standing tube)]	MS1760G	

Swabs shown under are in their actual size

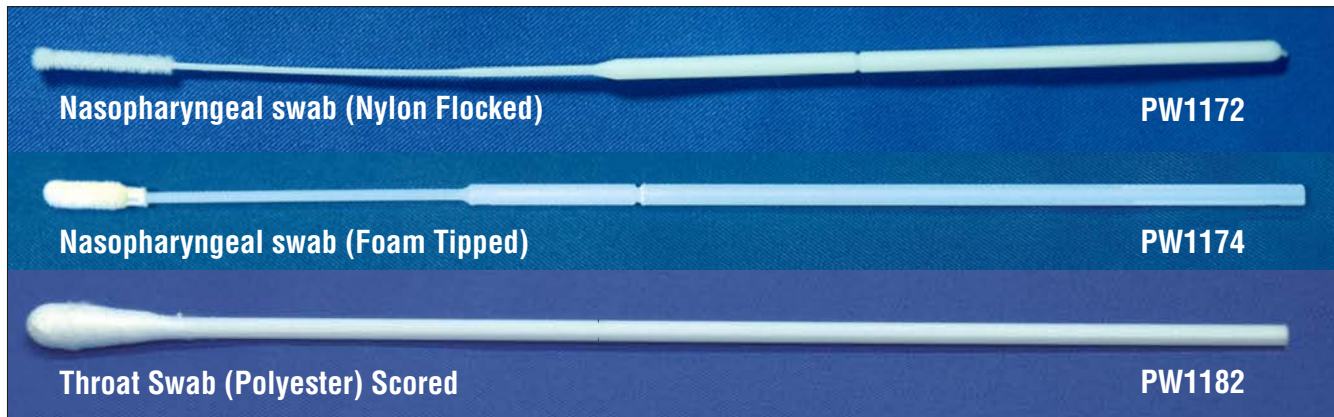

As Swine Flu Assumes epidemic proportions worldwide, and starts spreading its terror in India, HiMedia are ready with their plan to collect, transport and recover this virus under WHO approved guidelines for testing.

The only company worldwide to come out with the viable answer to the Mad Cow Disease scare with 1000 HiVeg Media, has now taken the lead to combat Swine Flu Viral Collection under international guidelines.



Production facility



R & D Centre



• • • expect only Quality from us™

HiMediaLaboratories™



HiMedia Laboratories Pvt. Limited

A-406, Bhaveshwar Plaza, Mumbai - 400 086, India.

Phone : 022-4095 1919 Fax : 022-4095 1920

Email : info@himedialabs.com

www.himedialabs.com

